SPARK Z4^a

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CONTENTS

SAFETY MEASURES	16
CLIMATIC CONDITIONS OF USE	16
MANUFACTURER'S WARRANTY	16
MAINTENANCE	16
OPERATING ELEMENTS	17
DESCRIPTION OF THE DEVICE	18
PACKAGE CONTENTS	28
TECHNICAL DATA	19
DISPLAY SYMBOLS	20
GETTING STARTED	20
OPERATION	21
1. Switching on the device	21
 The main menu	21
2. Switching off the device	21
3. Working with the generators of the device	21
 Submenu «Generators»	22
Sub-mode «Manual input»	22
 Sub-mode menu «Manual input»	22
Frequency sweep	23
Lock / Unlock / Pause	23
 Intermediate menu	23
Amplitude adjustment	23
Sub-mode «Compositions»	24
4. Programs on SD to a card	24
5. Internal memory	24
6. Settings	25
 Submenu «Settings»	25
Power supply	25
Warranty cards	26

SAFETY MEASURES

The device is made in accordance with the requirements for products with an internal power supply according to IEC 601-1-88.

Having discovered any problem when the device operates, immediately turn it off.

Do not use the device when its appearance is obviously damaged (broken case, broken screen).

Do not place the product near naked flames.

Do not leave the unit for long periods in direct sunlight.

Keep out of the reach of children.

Do not use the device in bathrooms and wet areas.

The application of the device for children is permitted only under adult supervision.

Do not place the device into a microwave oven.

Do not disassemble the device.

Do not immerse the device into water.

Do not drop the device.

Do not drop the other things on the device.

Protect the device from extreme temperature difference.

CLIMATIC CONDITIONS OF USE

Recommended ambient temperature from +10 ° C to +35 ° C.

The upper value of the relative air humidity 80% at +25 ° C.

Atmospheric pressure 620-800 mm Hg (84-106.7 MPa).

MANUFACTURER'S WARRANTY

1. The company - manufacturer guarantees that the device meets the requirements of Technical Specifications, provided the consumer meets the conditions of the use, transportation and storage.

2. The warranty period is 12 months from the date of putting into operation.


3. During the warranty period, the manufacturer replaces or repairs the broken-down device upon presentation of the guarantee card if there is no mechanical damage to the device.

4. The warranty does not cover the battery.

MAINTENANCE

Clean the external parts without the use of solvents and abrasives.

OPERATING ELEMENTS

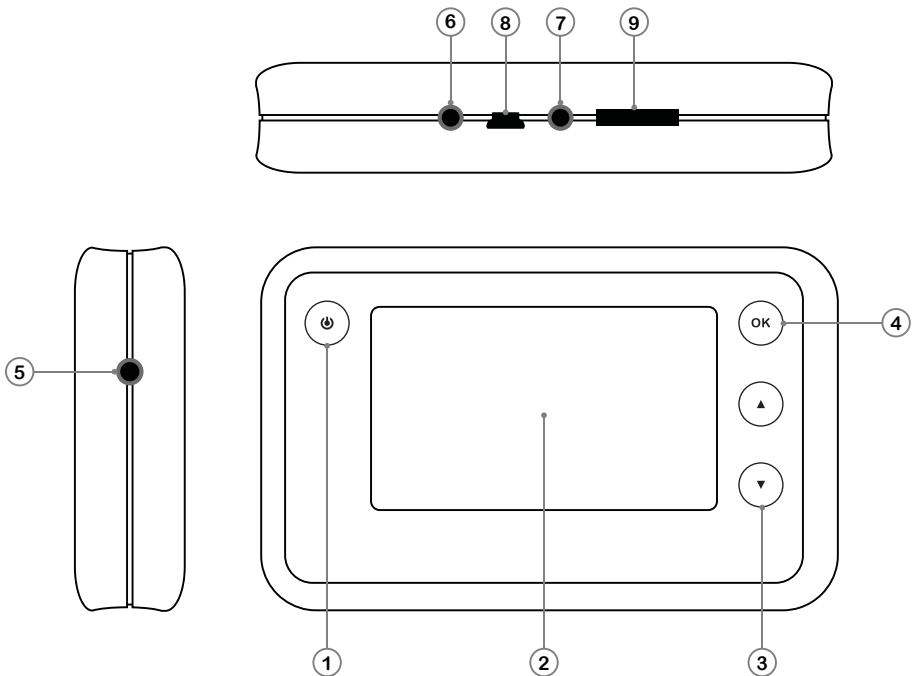- (1) On / off button of the device
- (2) Display
- (3) Menu navigation buttons
- (4) Choice confirmation button
- (5) Connector for the external inductor
- (6) Left electrode connector
- (7) Right electrode connector
- (8) Mini USB connector
- (9) Micro SD slot

The manufacturer reserves the right to change the configuration of the front panel, which will not involve the deterioration of the quality and usability of the device.

DESCRIPTION OF THE DEVICE

Device Spark Z4-a (hereinafter referred to as the device) is a portable electronic device programmed to generate low-intensity electromagnetic frequencies and induction field.

Two built-in independent generators:

- A generator with the output to the electrodes (range: 1 to 3,000,000 Hz), and the modulation option;
- A generator of low intensity electromagnetic fields with the output to the inductor (range: 1 to 25000 Hz);

It is possible to generate frequencies both with one of the generators, and simultaneously with two of them with independent frequency settings within the appropriate range.

For each generator there is a possibility to adjust the amplitude of pulse fluctuations and frequency "sweep" in the user-defined range.

The built-in memory of the device is used as alternate, in case of absence SD card which can be programmed by the user complexes. Otherwise, these complexes are written in, automatically generated to the SD card folder 'Personal.'

The support of MicroSD card up to 2GB significantly increases the number of recorded programs in the section «Programs on the SD card.»

Output frequencies for effects on the human body by using the pluggable electrodes and/or embedded in the housing coil. Optionally connect external loop-coil (sold separately). when connecting an external inductor alternator output automatically switches with built-in to.

The charging of the built-in rechargeable battery and connection to the computer are through the jack miniUSB.

PACKAGE CONTENTS

1	The device	1pc.
2	Cable	2pc
3	Cylindrical electrode	2pcs
4	Software CD (Optional)	1pc
5	5V AC adapter	1pc

TECHNICAL DATA

Current consumption, not more than	20 mA	
The range of generated frequencies	from 1 to 3,000,000 Hz	
Adjustment of the generators output amplitude	to 18 V	
Adjustment of the Electrodes generator signal intensity	9 levels	
Adjust. of the Inductor generator signal intensity	4 levels	
Modulation (for the electrodes generator)	to 100,000 Hz	
Frequency sweep range	from 1 to 3,000,000 Hz	
Display resolution	128 x 64 px	
Support for microSD memory cards up to 2GB	without the SDHC	
Dimensions	105x68x20 mm	
Net weight, not more than	0,112 kg	
Time of continuous operation of the device *	8 hours	
Internal battery	Capacity	1200 mAh
	Voltage	3.7 V
Charging the battery	Through USB	

* at a level of intensity not above the average or normal distribution of frequencies when the load is standard (50-100 ohms)

Warning!

If the device informs of the need to charge the built-in battery, it is necessary to charge the device, and in no case should you repeatedly try to turn it on at the same time. This would lead to failure of the battery itself.

When the generator operates at frequencies above 200 kHz, the current consumption increases and, accordingly, there is a reduction of the device operation time of one battery charge.

It is not advisable to use the electrical output of the generator for a long use with the resistance of the external load below 10 kilohms, for example, a container with liquid.

DISPLAY SYMBOLS

Symbol	Description
	Internal clock of the device (24-hour format)
	Generators' intensity indicator: the electrodes (E) and the inductor (I)
	Generators' open-circuit failure indication: the electrodes (E) and the inductor (I)
	The indication of generators' blocking by the user: the electrodes (E) and the inductor (I)
	Built-in inductor activation symbol
	The activation of the external inductor symbol (connected via connector 5)
	Only in the "Compositions". From the left to the right: program number and program cell number
	Program execution indication
	Group of programs execution indication
	Battery charge indicator
	Turned off loudspeaker indicator

GETTING STARTED

1. Pull out the device from the shipping container and with visual inspection make sure there is no mechanical damage.
2. After transportation or storage at low temperatures, the device must be kept at normal conditions for 2 hours.

OPERATION

1. Switching on the device.

Turn on the device by pressing and holding the button (1).

On the display (2) you see a greeting splash screen and after a few seconds the device switches to the interface language selection mode. A list of languages is offered. Navigate the menu with navigation buttons (3) and having chosen your language, click confirm button (4). Similarly, navigate the whole menu of the device.

The display shows the main menu, and voice prompt tells the device's readiness for operation.

Main Menu	Explanations
Generators	A sub-menu to set and run frequency generators
Built-in programs	A sub-menu to start programs installed by the manufacturer
Personal programs	A sub-menu to run programs stored via a USB port (8)
SD card programs	A sub-menu to run programs on a connected SD card (see paragraph 6)
Settings	A sub-menu to set the device (see paragraph 7)

2. Turning off the device.

Turn off the device by pressing and holding the button (1). If the device is in program run mode, you must first stop the program.

The device will turn off automatically if not to perform any action with any buttons within the time specified in the Settings -> Automatic turning off.

3. Working with generators of the device.

Get the device ready to work and enter the Generators Menu. The display shows the submenu "Generators".

Submenu "Generators"	Explanations
Manual input	Log in the sub-mode "Manual input"
Compositions	Log in sub-mode "Compositions". The programming of frequencies with a device
Back	Exit to the main menu (pressing the button (1) is always the exit to the upper level)

Sub-mode "Manual input"

In sub-mode Manual input two generators' independent settings are available, namely, Generator, with output to the Electrodes (the connectors (6) and (7), left and right, respectively), and the Inductor, with output to the external inductor or a loop (connected to the connector (5) located in the right-hand edge of the device).

Sub-mode menu "Manual input"		Explanations
Electrodes	On / Off	Enables / disables the output to the electrodes
The frequency	0000 000.00-3000 000.00 Hz	The frequency value of the Generator 1
Ranges to	0000 000.00-3000 000.00 Hz	The final value of the Generator 1
Modulation	000 000.00 - 100 000.00 Hz	Modulating frequency
Inductor	On / Off	On / Off access to the external inductor
The frequency	0000 000.00 - 25 000.00 Hz	The frequency value of Generator 2
Ranges to	0000 000.00 - 25 000.00 Hz	The final frequency value of Generator 2
Duration	00:01 - 59:59	The length of the program run (MM: SS)
Run		Running the generation program
Save as		Saving the parameters in the composition program cell
Back		Exit to the upper level

Navigating in the sub-mode menu, use the navigation buttons (3), using button (4) choose the line in which you want to set a parameter. After choosing the line, with the buttons (3) increase or decrease the value to the desired value and hold the button (4) to confirm and save the values. Also, a single press of the button (4) will move to a different number position.

Frequency sweep

A smooth frequency change occurs in the range between the first and second values. If the value "Change to" was not set, the frequency generation is constant. After setting the required generators' values, set the duration of the generation program and run the program. The program execution band and the remaining time will be shown on the display (2).

Lock / Unlock / Pause

During the program execution the buttons of the device are blocked to avoid accidental button pressing. To unlock the buttons, press buttons (1) and (4) at the same time.

If it is necessary to stop or pause the program, press the upper right and upper left buttons at the same time.

The program execution will pause. You will see an Intermediate menu on the display.

Intermediate menu	Explanations
Continue	Repair the open-circuit and continue the program
Turn off the electrodes or Turn off the inductor	If necessary, turn off the output
Cancel	Stopping the program

The intermediate menu is also displayed while there is an open-circuit. The program stops its work and shows a warning message.

Execute this necessary action.

At the end of the program, the display shows «Completed» and you can hear the voice message.

Double-clicking the button (4) resumes the program with the same settings. Single pressing of the button (1) or (4) returns you to the sub-mode menu «Manual input».

The tuned program can be stored in the composition program cell for the further usage from the device memory.

Amplitude adjustment

The amplitude of the output pulses to electrodes is adjusted by buttons (3).

The oscillation amplitude of the inductor output is adjusted with the same buttons, but with pre-pressed and held button (4).

The level of the controlled parameter is displayed in the top line of the display.

See display symbols.

The adjustments are possible only when the corresponding output is active and there is no open circuit of its load.

Sub-mode "Compositions"

Submode "Compositions" consists of a list of nine programs, each of these programs contains nine programmable cells. Cells are programmed with the help of a menu similar to the sub-mode menu "Manual input". The only difference between these menus - a line "Save". When you choose it, there is a rewriting of the current cell, if any changes have been made in the settings for the generators.

To create a composition, choose a program (one of nine), enter the "Cells" menu and adjust the parameters of the generators. If necessary, choose the next cell and adjust the new settings.

You can run the programmed cell from the current menu or the entire program, going to the upper level (the button (1) or "back").

In the list of programs, in front of each program you can see the total duration time of all cells recorded in this program. The program run is a consistent (according to the number) execution of cells frequencies.

4. Programs on SD to a card

In this section in a folder «Basic», the programs delivered with the device occur. And compiled user programs, with the help of specialized software, written to the SD card in the automatically generated 'personal' folder.

5. Internal memory

The internal memory of the device is used as a backup, to record Personal programs if the SD memory card is not being used by the user. Memory capacity approximately twenty Personal programs.

6. Settings

Adjust the device to suit your needs to use.

The Settings submenu	Explanatory Notes
Display	The settings of contrast, brightness, and the display backlight time
Sound	The setting of the audio signals and voice prompts volume (which can be switched off)
Automatic shut-off	The setting of the automatic shut-off interval when the device is idle
Time	The setting of the built-in device clock
Information	Displays the serial number of the device, firmware version and assembly
Back	Exit to the upper level (into the Main menu)

POWER SUPPLY

The device has a built-in rechargeable battery that requires periodic charging. In the case of low voltage below the allowed rate, the device will let you know with the help of the voice prompt, and if no measures are taken, it will turn off.

In the case of a strong battery discharge the device will not be able to turn on and then it is necessary to connect the charger or USB cable to the computer. The appropriate note will appear and the battery charging process will be shown on the display.

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